

Measurement of proton-driven antiport in Escherichia coli

Holdsworth, S. R., & Law, C. J. (2014). Measurement of proton-driven antiport in Escherichia coli. *Bio-protocol*, 4(21), [e1278]. <http://www.bio-protocol.org/e1278>

Published in:
Bio-protocol

Document Version:
Early version, also known as pre-print

Queen's University Belfast - Research Portal:
[Link to publication record in Queen's University Belfast Research Portal](#)

Publisher rights
Copyright 2015 The Author

General rights
Copyright for the publications made accessible via the Queen's University Belfast Research Portal is retained by the author(s) and / or other copyright owners and it is a condition of accessing these publications that users recognise and abide by the legal requirements associated with these rights.

Take down policy
The Research Portal is Queen's institutional repository that provides access to Queen's research output. Every effort has been made to ensure that content in the Research Portal does not infringe any person's rights, or applicable UK laws. If you discover content in the Research Portal that you believe breaches copyright or violates any law, please contact openaccess@qub.ac.uk.

Category by field: Microbiology > Microbial biochemistry > Other compound
Category by field: Microbiology > Microbial metabolism > Nutrient transport
Category by organism: Bacteria > *Escherichia* > *Escherichia coli* > Physiology

Measurement of Proton-driven Antiport in *Escherichia coli*

Scarlett R. Holdsworth¹ and Christopher J. Law^{1*}

¹Institute for Global Food Security, School of Biological Sciences, Queen's University Belfast, 97
Lisburn Road, Belfast BT9 7BL, United Kingdom

*Corresponding author: Christopher J. Law, c.law@qub.ac.uk

[Abstract] Secondary active transport of substrates across the inner membrane is vital to the bacterial cell. Of the secondary active transporter families, the ubiquitous major facilitator superfamily (MFS) is the largest and most functionally diverse (Reddy *et al.*, 2012). Recently, it was reported that the MFS multidrug efflux protein MdtM from *Escherichia coli* functions physiologically in protection of bacterial cells against bile salts (Paul *et al.*, 2014). The MdtM transporter imparts bile salt resistance to the bacterial cell by coupling the exchange of external protons (H⁺) to the efflux of bile salts from the cell interior via an antiport reaction. This protocol describes, using fluorometry, how to detect the bile salt/H⁺ antiport activity of MdtM in inverted membrane vesicles of an antiporter-deficient strain of *E. coli* TO114 cells by measuring transmembrane Δ pH. This method exploits the changes that occur in the intensity of the fluorescence signal (quenching and dequenching) of the pH-sensitive dye acridine orange in response to changes in [H⁺] in the vesicular lumen. Due to low levels of endogenous transporter expression that would normally make the contribution of individual transporters such as MdtM to proton-driven antiport difficult to detect, the method typically necessitates that the transporter of interest be overexpressed from a multicopy plasmid. Although the protocol described here pertains to measurement of bile salt efflux by a specific membrane transporter protein (MdtM), it can be easily adapted for measurement of antiport of other substrates by any other antiporter that exchanges protons for countersubstrate. We have therefore omitted description of the bacterial growth and protein overexpression conditions (as these are usually different for every membrane protein), and begin the procedure from the point of harvesting of 1 L cultures of bacterial cells enriched with plasmid-encoded antiporter.

Materials and Reagents

1. *Escherichia coli* (*E. coli*) TO114 cells (gift of Prof. Hiroshi Kobayashi, Chiba University, Japan)
2. Acridine orange hemi (zinc chloride) salt (Sigma-Aldrich, catalog number: A6014)
3. BisTris propane (BTP) (Sigma-Aldrich, catalog number: B6755)
4. Phenylmethanesulfonyl fluoride (PMSF) (Sigma-Aldrich, catalog number: P7626)
5. Deoxyribonuclease I (DNase) from bovine pancreas (Sigma-Aldrich, catalog number: DN25)
6. Carbonyl cyanide 3-chlorophenylhydrazone (CCCP) (Sigma-Aldrich, catalog number: C2759)

7. Sodium DL-lactate solution 50% aqueous (VWR International, catalog number: 27927.298)
8. Magnesium sulphate heptahydrate ($\text{MgSO}_4 \cdot 7\text{H}_2\text{O}$) (Thermo Fisher Scientific, catalog number: M/1000/60)
9. Sodium cholate hydrate (Sigma-Aldrich, catalog number: C1254)
10. Absolute ethanol (Thermo Fisher Scientific, catalog number: E/0650DF/17)
11. High purity (18 M Ω) Millipore or AnalR water
12. Choline chloride (Sigma-Aldrich, catalog number: C1879)
13. Sucrose (Sigma-Aldrich, catalog number: 84097)
14. DL-dithiothreitol (Sigma-Aldrich, catalog number: 43815)
15. TRIZMA base (Sigma-Aldrich, catalog number: T1503)
16. 32% hydrochloric acid (VWR International, catalog number: 20254.321)
17. Tris/choline/dithiothreitol/sucrose (TCDS) buffer (see Recipes)
18. Transport assay buffer (see Recipes)
19. 200 mM sodium DL-lactate solution (see Recipes)

Equipment

1. Temperature-controlled shaking incubator for bacterial growth
2. Large ice bucket
3. Refrigerated, large capacity centrifuge and rotor for harvesting bacterial cells
4. 1,000 ml or 500 ml centrifuge pots and lids
5. Benchtop vortexer
6. 100 ml conical flask
7. 100 ml beaker and stir bar to fit
8. Magnetic stirrer
9. 25 ml disposable plastic pipettes
10. Selection of single channel pipettes (1,000 μl , 100 μl , 20 μl , 10 μl)
11. Pipette tips for above
12. Refrigerated centrifuge and rotor capable of spinning ~50 ml tubes at 18k x g
13. Refrigerated ultracentrifuge, rotor and polycarbonate ultracentrifuge tubes capable of handling ~30 - 50 ml volumes
14. French Press (Thermo Fisher Scientific, catalogue number: FA-078)
15. Standard pressure cell (40 kpsi; 35 ml capacity) (Thermo Fisher Scientific, catalog number: FA-031)
16. 50 ml syringes (for filtering solutions)
17. 0.22 μm sterile filters to fit 50 ml syringe
18. 1 L and 500 ml Duran bottles with lids
19. 1.5 ml Eppendorf tubes
20. Medical wipes (Kimwipes®)

21. Parafilm
22. UV/vis spectrophotometer and 10 mm pathlength quartz cuvette
23. 10 x 4 mm, 1,400 µl volume quartz cuvette for fluorescence spectroscopy (Hellma, catalog number: 104F-QG)
24. Small magnetic stir bar to fit inside quartz cuvette
25. Fluorometer *e.g.* Fluoromax-4 (Horiba) capable of performing time-based acquisition measurements and fitted with a temperature controlled cuvette holder and stirrer
26. Source of compressed air (either aerosol can or fixed supply) for drying cuvette

Procedure

A. Harvesting bacterial cultures and preparation of inverted membrane vesicles

1. Transfer the *E. coli* TO114 cells that contain overexpressed transporter into 500 ml or 1,000 ml capacity centrifuge pots that have been pre-chilled on ice for at least 15 mins. Pre-cool the centrifuge to 4 °C and harvest cells by centrifugation at 5,000 x *g* for 20 min.
2. Decant the supernatant and wash the pelleted cells by resuspending in chilled (4 °C) TCDS buffer. Use 30 ml of TCDS buffer for each one litre of cell culture that was pelleted. Maintain the cells on ice during this procedure. Harvest the washed cells by centrifugation as described in step A1 (above) and repeat the washing procedure. Resuspend the resultant cell pellet in 30 ml of chilled TCDS buffer containing 2 mM PMSF (Which should be made up as a 100 mM stock in ethanol and stored at -20 °C; ensure the solution is thawed thoroughly and vortexed vigorously before use.) and 5 µM DNase and maintain the mixture on ice.

Notes:

- a. Typically, a 1 L culture of bacterial cells will provide sufficient material for these experiments.*
- b. Cells should be resuspended either by gentle vortexing using a benchtop vortexer or by gentle aspiration using a 25 ml sterile plastic pipette.*
3. Decant the resuspended cells into a 100 ml beaker containing a stir bar, place on magnetic stirrer and stir in a cold room at 4 °C or on ice for 20 min.
4. Generate inverted membrane vesicles by a single passage of the resuspended cell mixture through a French pressure cell at a minimum of 4,000 psi. If the pressure is too low, inverted vesicles will not be formed. The pressure cell should be chilled on ice for ~30 min prior to use. The vesicle mixture should be collected in a 100 ml conical flask kept on ice.
5. Decant the mixture into a pre-chilled ~50 ml centrifuge tube and remove any unbroken cells and cell debris by centrifugation at 18,000 x *g* for 10 min at 4 °C. Carefully decant the supernatant containing the cell membrane vesicles into a pre-chilled 30-50 ml volume polycarbonate ultracentrifuge tube on ice.

6. Harvest the inverted vesicles by ultracentrifugation at 100,000 x g for 1 h at 4 °C. Carefully decant the supernatant and retain the pellet. Place the ultracentrifuge tube containing the pelleted vesicles on ice.
7. Resuspend the inverted vesicle pellet in 1 ml of ice-cold TCDS buffer by gentle aspiration using a 1,000 µl pipette. Transfer the resuspended vesicles to a pre-chilled 1.5 ml Eppendorf tube on ice for use in the transport assay. In our experience, vesicles stored on ice are stable for several hours.
8. Quantify the total membrane protein content of the inverted vesicles by UV absorbance spectroscopy at 280 nm. Blank the spectrophotometer using a 10 mm pathlength quartz cuvette containing 1,000 µl of TCDS buffer. The buffer should be at room temperature to prevent frosting of the cuvette faces. Clean the faces of the cuvette using a fresh paper wipe prior to measurement. Once the spectrophotometer is blanked, remove 5 µl of buffer from the cuvette using a 10 µl pipette and replace with 5 µl of vesicles. Cover the opening of the cuvette with a square of Parafilm and invert the cuvette a few times to ensure a homogeneous distribution of vesicles. Record the absorbance of the vesicle mixture at 280 nm and calculate the total membrane protein concentration assuming that an A_{280} of 1.0 is equivalent to a protein concentration of 1.0 mg/ml.

Notes:

- a. Remember to multiply the 280 nm absorbance value you obtain by a factor of 200 to calculate the concentration of the undiluted vesicle mixture in mg/ml.
- b. At this stage of the preparation, the resuspended vesicles can be transferred to tubes in aliquots of 25-100 µl, snap-frozen in liquid nitrogen and stored at -80 °C for subsequent use. Vesicles frozen in this way will retain their integrity for several months. However, if frozen vesicle stocks are used for the subsequent transport measurements, the vesicles must be thawed very slowly on ice prior to use to prevent their fracture.

B. Fluorometric antiport assays

1. In this section we describe the set-up parameters for a Fluoromax-4 fluorometer. These parameters, however, can form the basis for the set-up of other fluorometers. Once the instrument is switched on and the software booted up, set the temperature of the cuvette holder to 25 °C. Open the instrument software and select for time-based data acquisition with excitation and emission wavelengths of 492 nm and 525 nm, respectively. Set the excitation and emission slit widths to 1.5 nm and 2.5 nm, respectively.
2. Add an aliquot of inverted vesicles (which should be maintained on ice in TCDS buffer) to room temperature transport assay buffer containing the acridine orange probe in a 10 mm x 4 mm quartz cuvette to a final concentration of 0.5 mg/ml membrane protein in a total volume of 1,500 µl. The longest pathlength of the cuvette should face the excitation light source. Place a small magnetic flea into the cuvette and stir the contents gently. Allow the vesicles and assay buffer to equilibrate for ~200 sec.

3. Start recording the fluorescence emission. After approximately 50 sec, add 15 μ l of 200 mM stock sodium DL-lactate solution to the cuvette contents to give a final sodium DL-lactate concentration of \sim 2.0 mM. Addition of lactate energises the vesicles and generates a respiration dependent Δ pH (acid inside) across the inverted vesicle membrane as H^+ is pumped into the vesicle interior. This causes a dequench of the acridine orange fluorescence signal (see Figure 1 in Representative data).
4. Following the establishment of a Δ pH, monitor the acridine orange fluorescence dequench for a further \sim 200 sec until it stabilises. Initiate MdtM-mediated, proton-driven antiport by adding substrate (in this case the bile salt sodium cholate) to the inverted vesicle mixture. We added 12.5 μ l of 250 mM stock solution made up in high purity water to give a final concentration of sodium cholate in the cuvette of \sim 2.0 mM. For other substrates we suggest testing a range varying from 1 mM to 100 mM to establish the concentration that gives the best dequench signal. Upon addition of substrate, there should be an immediate dequench (rise) of the acridine orange fluorescence emission signal (see Figure 1a) due to dissipation of the established Δ pH as a result of MdtM-mediated sodium cholate/ H^+ antiport activity and concomitant alkalinisation of the vesicle lumen.
5. Record the fluorescence dequench signal for 50 sec to allow the antiport reaction to achieve a steady state (as observed by a plateauing of the fluorescence signal). Addition of the protonophore CCCP to a final concentration of \sim 100 μ M (1.6 μ l of a 100 mM stock made up in ethanol) in the assay mixture should be performed to abolish Δ pH-driven, MdtM-mediated antiport activity. Addition of CCCP will cause a further dequench of the fluorescence signal as the Δ pH is dissipated (see Figure 1). Record the fluorescence signal for a further \sim 50 sec then terminate the acquisition and save the electronic data.
6. Decant the cuvette contents into a suitable waste container and wash the cuvette thoroughly with ethanol then high-purity water. Dry the cuvette carefully with compressed air.

Representative data

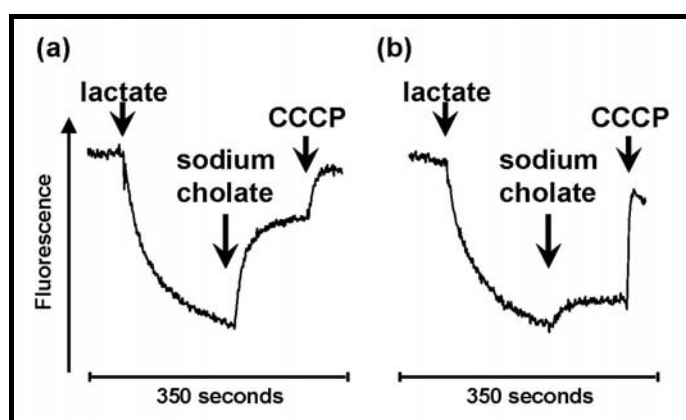


Figure 1. Representative measurements of the fluorescence quench/dequench of acridine orange upon addition of bile salts to inverted vesicles of *E. coli* TO114 cells that

overproduced recombinant (a) wild type MdtM or, as a control, (b) the dysfunctional MdtM D22A mutant. Respiration-dependent generation of Δ pH (acid inside) was established by addition of sodium DL-lactate as indicated. Sodium cholate was added to vesicles as indicated to initiate the transport reaction and CCCP was used to dissipate Δ pH. The fluorescence dequench observed in the control experiment (panel b) upon addition of sodium cholate is due to antiport activity of chromosomally encoded MdtM. The traces are representative of experiments performed in triplicate on at least two separate preparations of inverted vesicles.

Notes

- a. *To ensure reproducibility, the assays should be performed in triplicate on at least two separate preparations of inverted vesicles.*
- b. *If the fluorescence signal does not quench, or enhances, the vesicles have either not maintained integrity or are not inverted and their preparation needs to be repeated.*
- c. *As with all assays that rely on detection of fluorescence, robust controls must be in place to ensure that any detected transport activity can be attributed unambiguously to the protein of interest. In our experiments, we used inverted vesicles that overexpressed MdtM D22A, a dysfunctional point mutant of MdtM, as a negative control (see Figure 1b).*
- d. *If the method is to be used for detection of metal ion/ H^+ antiport activity, the use of inverted vesicles generated from the antiporter-deficient TO114 strain of *E. coli* is important because at least four other transporters (NhaA, NhaB, ChaA and MdfA) present in the bacterium catalyse a monovalent metal cation/ H^+ exchange.*
- e. *Finally, if this protocol is used for comparison of antiport activities of wild type and mutant transporters, the amount of target protein present in the inverted vesicle membranes must be quantified (usually by immunodetection methods) to ensure that any measured differences in H^+ uptake are due solely to differences in the activity of the transporters and not to differences in expression levels.*

Recipes

1. Tris/choline/dithiothreitol/sucrose (TCDS) buffer (1 L)
Consisting of 10 mM Tris-HCl (pH 7.5), 140 mM choline chloride, 0.5 mM dithiothreitol and 250 mM sucrose
1 M Tris-HCl (pH 7.5) 10 ml
2 M choline chloride 70 ml
1 M sucrose 250 ml
Make up to 999.5 ml with high purity water then check pH and adjust if necessary with HCl
Autoclave or sterile filter and store at 4 °C
Add 0.5 ml dithiothreitol from frozen 1 M stocks immediately before use

2. Transport assay buffer (100 ml)
Consisting of 10 mM BisTris propane (pH 7.2), 5 mM MgSO₄, 1 μM acridine orange
100 mM BisTris propane (pH 7.2) 10 ml
2 M MgSO₄ 0.25 ml
10 mM acridine orange 10 μl
Make up to 100 ml with high purity water then sterile filter
Acridine orange is light sensitive so both the transport assay buffer and acridine orange
stocks should be stored either in amber bottles or in a container protected from light at 4 °C in
the dark.
3. 200 mM sodium DL-lactate solution (10 ml)
0.45 ml 50% sodium DL-lactate aqueous solution
Make up to 10 ml with high-purity water
Sterile filter and stored at 4 °C

Acknowledgements

This work was supported by BBSRC grant BB/K014226/1 (to CJL). The protocol described above
is adapted from one reported previously (Resch *et al.*, 2010).

References

1. Reddy, V. S., Shlykov, M. A., Castillo, R., Sun, E. I. and Saier, M. H., Jr. (2012). [The major facilitator superfamily \(MFS\) revisited](#). *FEBS J* 279(11): 2022-2035.
2. Paul, S., Alegre, K. O., Holdsworth, S. R., Rice, M., Brown, J. A., McVeigh, P., Kelly, S. M. and Law, C. J. (2014). [A single-component multidrug transporter of the major facilitator superfamily is part of a network that protects *Escherichia coli* from bile salt stress](#). *Mol Microbiol* 92(4): 872-884.
3. Resch, C. T., Winogrodzki, J. L., Patterson, C. T., Lind, E. J., Quinn, M. J., Dibrov, P. and Hase, C. C. (2010). [The putative Na⁺/H⁺ antiporter of *Vibrio cholerae*, Vc-NhaP2, mediates the specific K⁺/H⁺ exchange *in vivo*](#). *Biochemistry* 49(11): 2520-2528.